



High dimensional biology with high throughput multiomic profiling of single cells

Aruna Ayer¹, Hye-Won Song², Larry Wang¹, Devon Jensen¹, Chelsea Gordon¹, Yan Chen¹, Jamie Moskwa¹, Jody Martin¹, Tracy Campbell¹, Adam Wright², Zhiyuan Han¹, Gracia Genero¹, Raghavendra Padmanaban¹, Hongdian Huang¹, Quyen Bao¹, Ruifang Li¹, Thomas McCarthy¹, Rosary Nguyen¹, Ricelle Acob¹, Cynthia Sakofsky¹, Manish Takran¹, Samatha Vadreru¹, Joe Olives², Elham Hatami², Zhiqi Zhang², Samuel Shi¹, Punya Narayan¹, Lin Chen¹, Anne Tran¹, Yingwah Lee¹, Amie Radenbaugh¹, Youngsook Kim¹


¹BD Biosciences 2350 Qume Drive, San Jose CA; ²BD Biosciences 10975 Torreyana Rd, San Diego CA

Abstract


Flow cytometry technology remains one of the most powerful tools used by life scientists to analyze single cells and characterize the cell at the protein layer. As biologists have sought to know more, the single "ome" of proteins alone could not fully explain cellular complexity. Single cell multiomics provides a comprehensive view of the cellular heterogeneity and the complex interplay between multiple layers of cellular "omics", namely the DNA, RNA and protein layers. Researchers are increasingly applying single cell multiomics across many fields including cancer research, drug discovery, infectious disease research and more.

Here in, we present data to showcase how BD Rhapsody HT-Xpress System fits within a reliable end-to-end single-cell multiomic solution with upstream cell-sorting solutions and downstream single-cell multiomic assay kits and bioinformatic tools. The gentle and robust micro-well based high throughput single-cell partitioning system featuring a flexible 8-lane cartridge design that can process up to half a millions cells per cartridge. Cell capture rates are typically > 80%, even with cells of varying sizes and fragility, including neutrophils and natural killer cells and even nuclei. The comprehensive single-cell multiomic assays include DNA chromatin accessibility assay (ATAC-seq), targeted/whole transcriptome assay (WTA), TCR/BCR full-length sequencing assay and CITE-seq assay with Abseq antibody oligos for surface and intracellular protein profiling. We present data highlighting the sensitivity and specificity of these multiomic assays and showcase the depth of molecular information and discovery enabled by single-cell multiomics.

End to end single-cell multiomic workflow

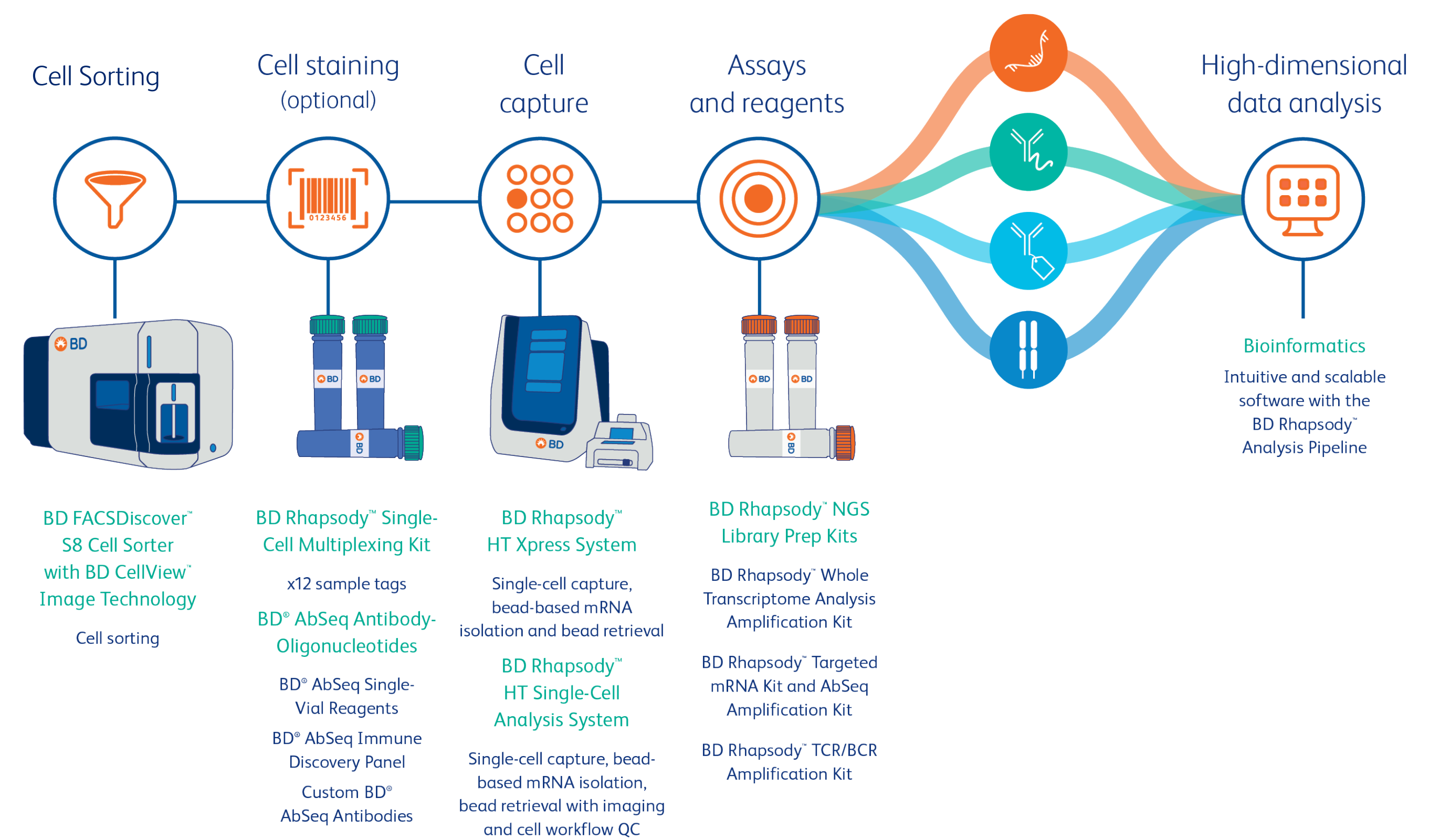

BD Rhapsody™ - Flexible, high-throughput platform for single-cell multiomics

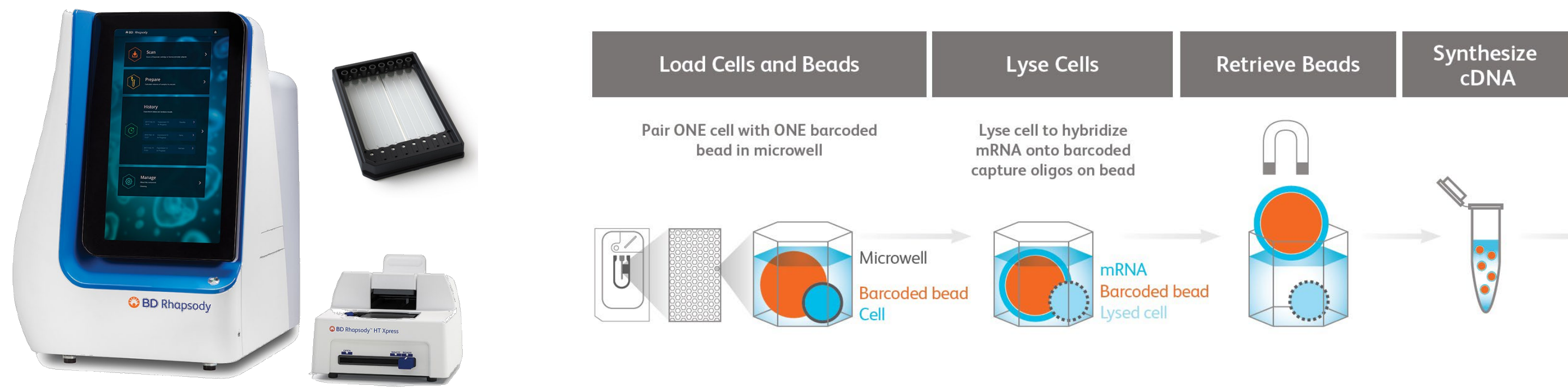


- Flexible, higher throughput cell capture with gentle micro-well technology
 - Load up to 500,000 cells per 8-lane cartridge, 100-65K cells per lane.
 - Partial use cartridge stable up to 6 months
 - cDNA on beads stable up to 1 year.
- Whole/Targeted transcriptome (mRNA), CITE-seq (mRNA+surface/intracellular protein), TCR/BCR Multiomic Assay (VDJ+mRNA+surface protein) assays, ATACseq (Standalone or multiomic chromatin DNA/mRNA Assay)
- No batch effect between lanes and cartridges.
- No sample loss due to clogging

BD Rhapsody™ TCR/BCR Next Multiomic Assay


Multiomic TCR/BCR sequencing assay with whole transcriptome and surface protein profiling

Next assay launching soon with improved full-length TCR pairing efficiency compared to V1.0 assay

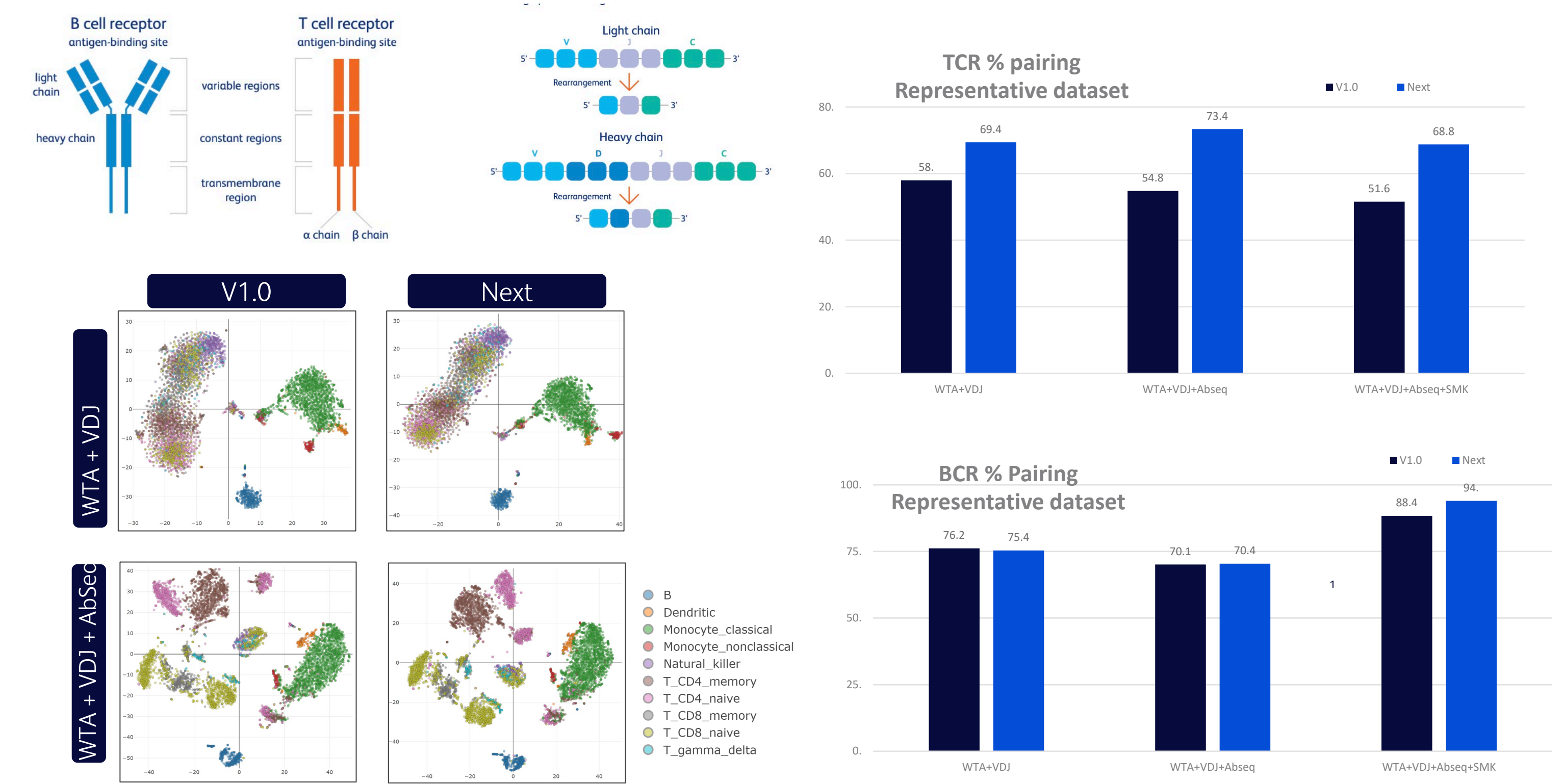


BD Rhapsody™ ATAC-Seq Assay

Standalone or multiomic single-nuclear chromatin accessibility assay with whole transcriptome

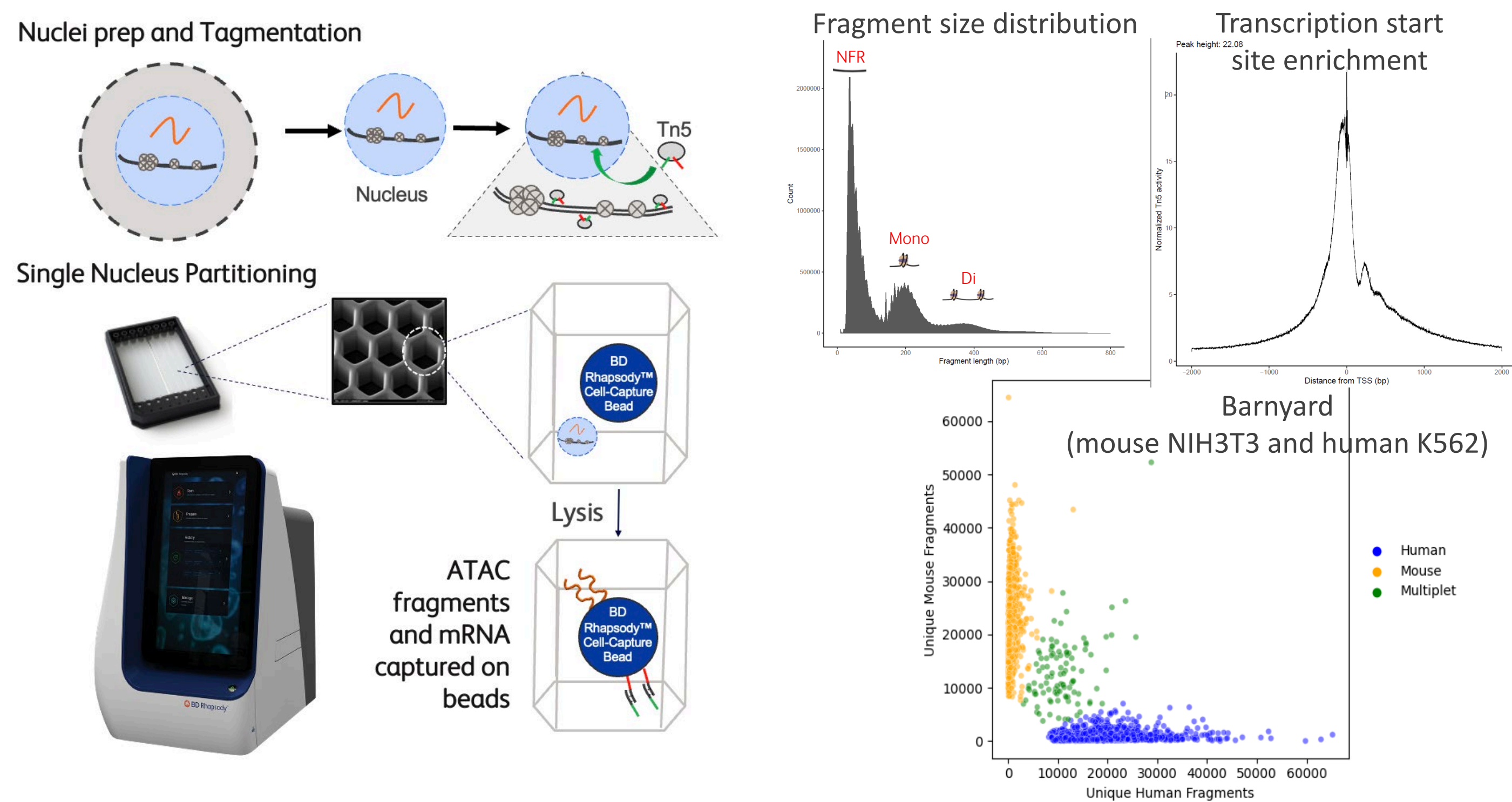

High specificity and sensitivity assay to profile open chromatin regions of the genome

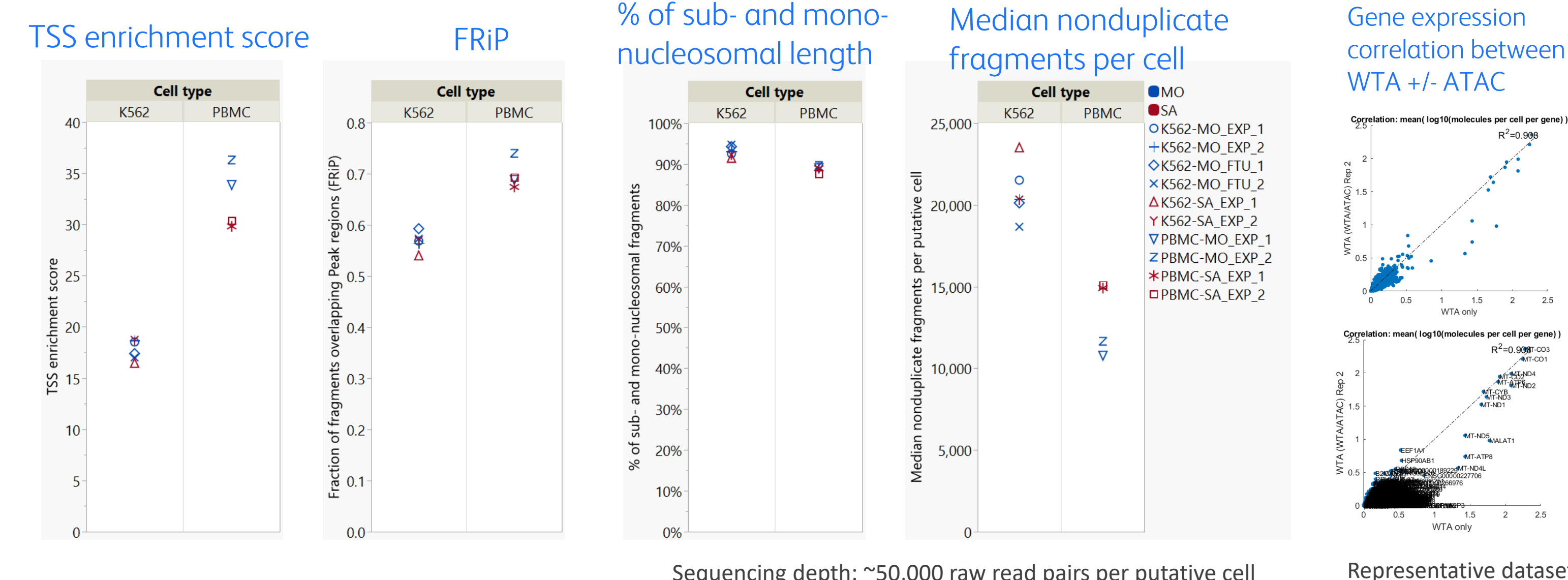


Conclusions

In summary, we show case how high parameter spectral flow cytometry with real-time imaging and flexible/high-throughput single-cell multiomic assays to profile DNA/RNA/protein enable deeper biological insights.

BD Rhapsody™ Intracellular CITE-seq Assay

Simultaneous profiling of RNA with surface and intracellular proteins on the BD Rhapsody™ System

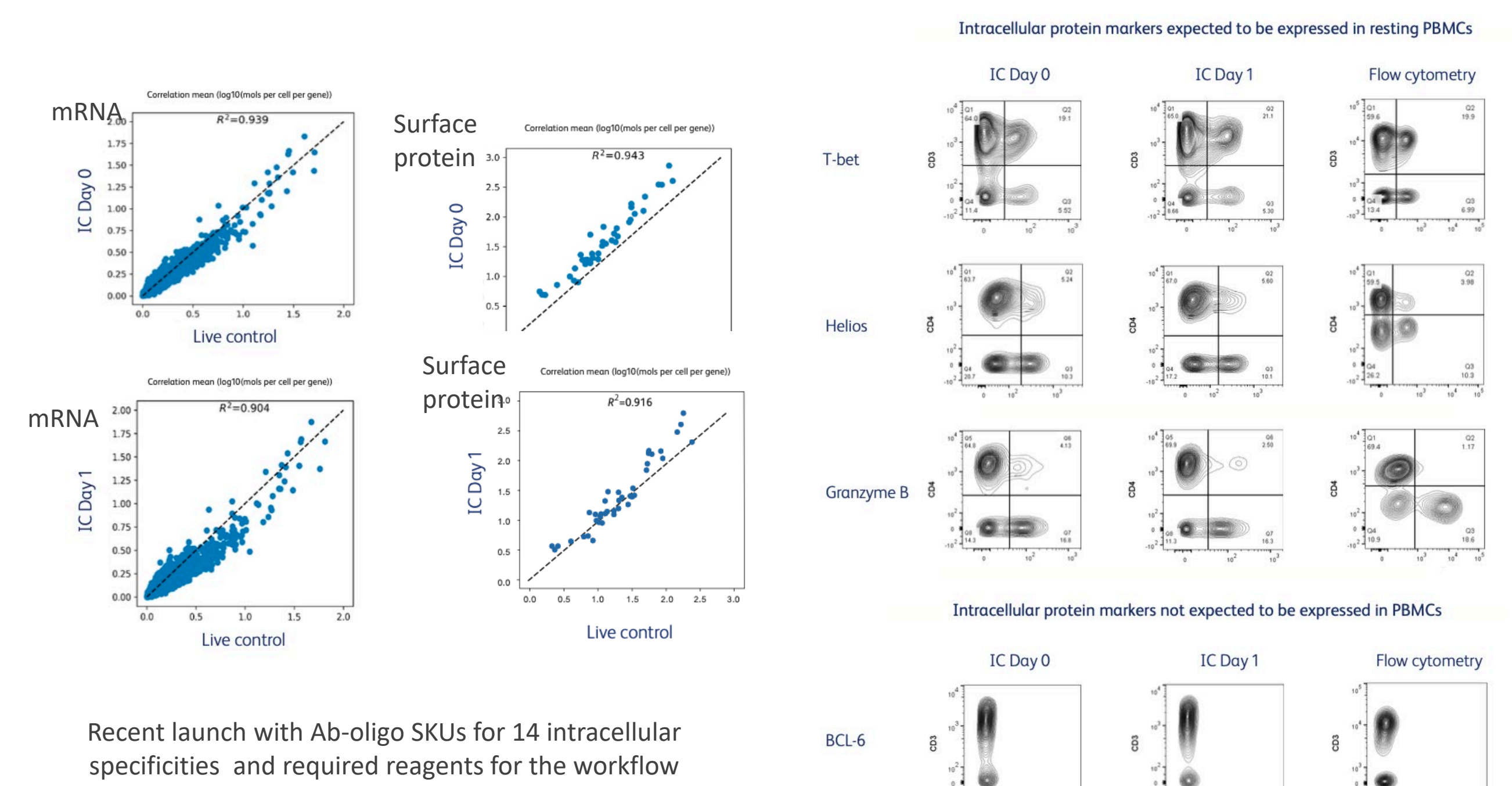


BD OMICS-Guard™ Sample Preservation buffer

PFA-free buffer that preserves RNA and protein epitopes in PBMCs for up to 72 h at 4 °C

- Biological sample preservation is a critical need
- De-centralized sample collection
- Longer sample processing times
- Large-scale population studies

- Available in two formats:
 - 50-mL bottle
 - 12 x 1-mL vials

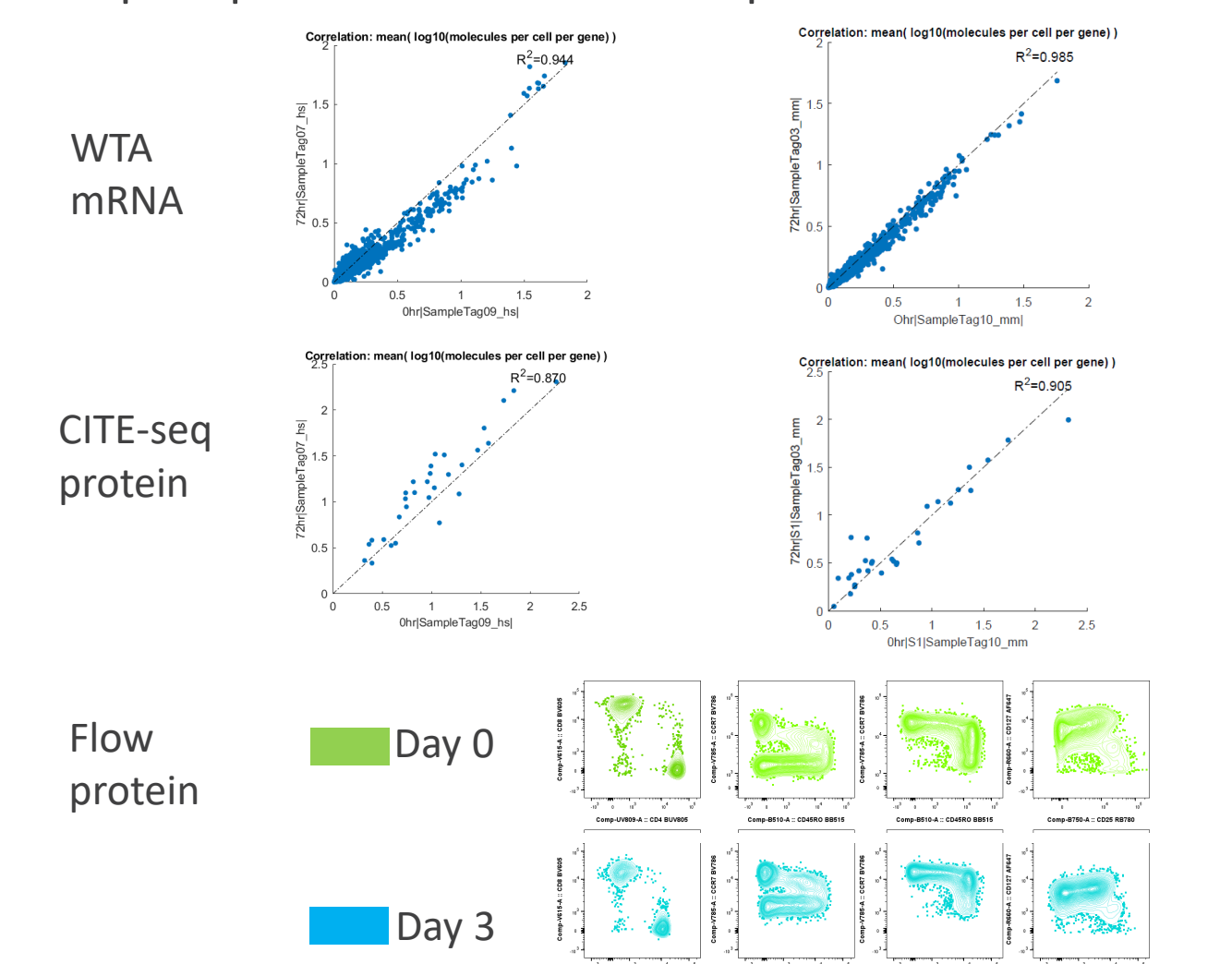

Class 1 Laser Product. For Research Use Only. Not for use in diagnostic or therapeutic procedures.

BD and the BD logo are trademarks of Becton, Dickinson and Company or its affiliates. All other trademarks are the respective property of their owners. © 2024 BD. All rights reserved. BD-120943(v1.0) 0524